



# PRODUCT CATALOG

THROW AWAY  
CONTRACTS  
FOREVER

DRAMATICALLY  
IMPROVE  
YOUR BOTTOM LINE

NEVER  
COMPROMISE  
ON QUALITY

RELIABLE  
SUPPLY  
PARTNERSHIP



YOUR PARTNER  
FOR GI ENDOSCOPY  
DEVICES AND ACCESSORIES



# TABLE OF CONTENTS

BIOPSY	6
DILATION	11
HEMOSTASIS	18
INFECTION PREVENTION	26
POLYPECTOMY	30
GI ACCESSORIES	37
PRODUCT REFERENCE CHART	41

## MICRO-TECH ENDOSCOPY IS PROUD TO OFFER YOU A COMPLETELY UNIQUE, COMPLETELY EXCEPTIONAL WAY OF DOING BUSINESS.

Imagine working with an organization that puts all its resources into focusing on you. A company that gives you the freedom to choose on your own terms, and will never lock you into a contract.

We operate by a higher standard. We know that you're busy, value efficiency, and require perfection from the devices you use to keep your patients healthy.

At Micro-Tech, you have the freedom to purchase what you need, when you need it. The highest quality will never be compromised, and your satisfaction is always guaranteed.

With facilities in Asia, Europe, and America, we leverage our global reach to continually refine and provide the best products for endoscopic care.

**You provide exceptional care. We provide exceptional support. Why settle for anything less?**



Micro-Tech Endoscopy has teamed up with Henry Schein Medical in an exclusive, innovative partnership to deliver more value to our customers. You can now combine your Micro-Tech Endoscopy and Henry Schein orders for a true one-stop shopping experience.

CURRENTLY SERVING OVER  
**3,500 ASCs**

NEXT DAY DELIVERY\*  
LATE DAY ORDERING CUT-OFF TIMES  
NOMINAL SHIPPING AND HANDLING CHARGES\*\*  
ORDER ITEMS IN LOW UNITS OF MEASURE†

## HOW TO ORDER



### CUSTOMER SERVICE

1.800.772.4346  
8:00 AM – 8:30 PM EST



### EMAIL

letsmeet@henryschein.com



### LOGIN OR CREATE AN ACCOUNT

www.henryschein.com



### FOR PRODUCT SUPPORT

1.800.447.1320



### EMAIL

info@mymicro-tech.com

\*For majority of customers – some restrictions apply

\*\*Minimum purchase of \$200.00

†Excludes bite blocks





# BIOPSY

LARGE CAPACITY BIOPSY FORCEPS

GI CYTOLOGY BRUSHES

BIOPSY FORCEPS

RETRIEVAL FORCEPS



## LARGE CAPACITY BIOPSY FORCEPS

Twin-wire mechanism and precision stamping technology for assured tissue sampling

- Pivoting jaws to aid in sampling where space may be limited
- Elongated jaw for greater tissue retention
- With needle option for multiples bites in a single pass
- Entry - Exit indicator markings on sheath
- Flexible and pushable smooth catheter

### BIOPSY FORCEP FOR 2.8MM SCOPE CHANNELS STERILE

Order Number	Henry Schein Item Number	Description	Jaw Style	Sheath Color	Working Length (cm)	Minimum Channel Size (mm)	Package Units
BF40286	133-4667	Large Capacity Forceps	Alligator	Blue	230	2.8	50/Box
BF40306	133-4666	Large Capacity Forceps with Spike	Alligator	Blue	230	2.8	50/Box



## GI CYTOLOGY BRUSHES

Stiff brushes for reliable tissue sampling

### CYTOLOGY BRUSH STERILE


Order Number	Henry Schein Item Number	Description	Brush O.D. (mm)	Catheter O.D. (mm)	Working Length (cm)	Minimum Channel Size (mm)	Package Units
CY49021	133-5541	Cytology Brush	2.0	1.8	180	2.0	20/Box
CY49051	133-5542	Cytology Brush	3.0	2.3	240	2.8	20/Box






# BIOPSY FORCEPS

With our wide selection of cost-effective disposable biopsy forceps, you are perfectly equipped for every procedure.

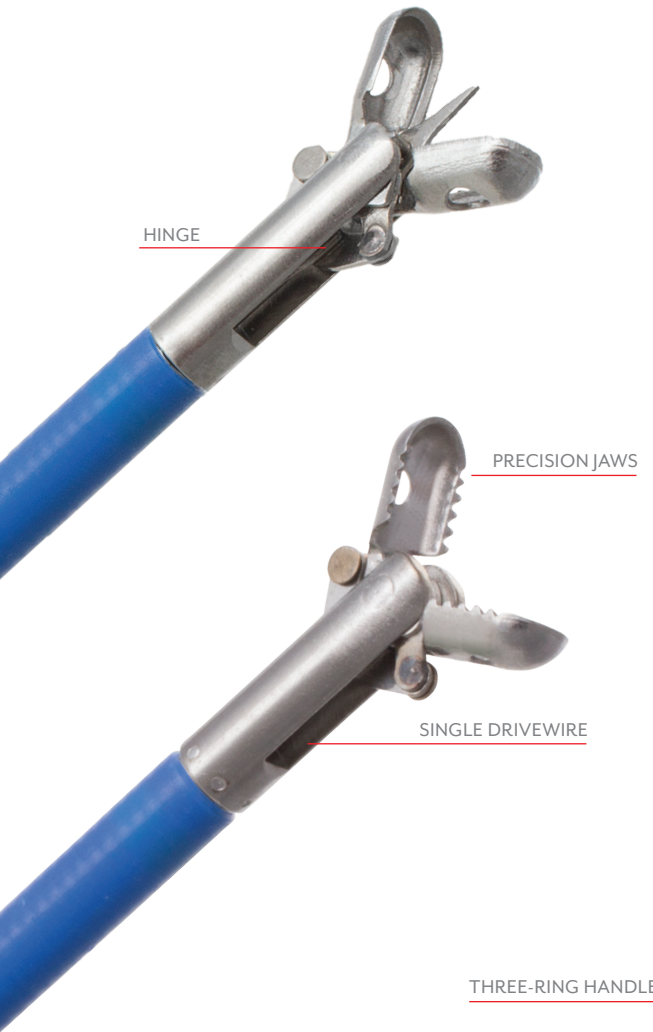

Available jaw diameters of 1.8 mm, 2.3 mm and 3.0 mm, and are offered in Oval, Serrated, Alligator jaws, with or without a spike  
 Available lengths of 160 cm, 180 cm and 230 mm  
 Standard Oval and Serrated offer larger capacity  
 All forceps feature entry and exit indicator markings

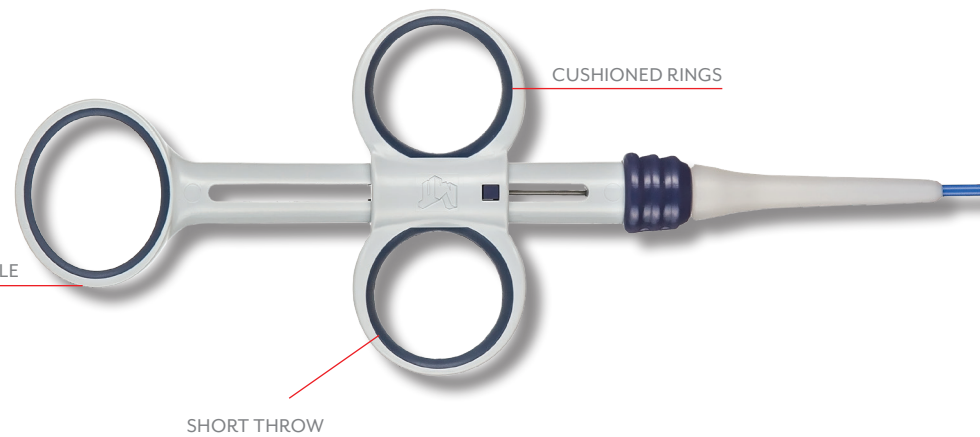

## KEY BENEFITS

### ERGONOMIC HANDLE

Three-ring comfort with cushioned lining and short throw provides improved grip and accommodates smaller hands.

### DIRECT DRIVE SYSTEM

The single drive wire and hinge linkage delivers smooth activation around curves and tortuous anatomy. The system delivers more closure force to the jaws which provides for improved cutting ability and assures complete jaw closure.

### SCOPE FRIENDLY DESIGN

The design of the Micro-Tech forceps family is based upon the robust and proven design of original reusable biopsy instruments. These designs provide optimal bite, while minimizing the risk of scope damage.

## BIOPSY FORCEP FOR 2.0MM SCOPE CHANNELS STERILE ②

Order Number	Henry Schein Item Number	Description	Jaw Style	Sheath Color	Working Length (cm)	Minimum Channel Size (mm)	Package Units
BF40021	129-0175	Mini Forceps	Oval	Yellow	160	2	10/Box
BF40051	129-0176	Mini Forceps with Spike	Oval	Yellow	160	2	10/Box
BF40141	129-0181	Mini Forceps	Alligator	Yellow	160	2	10/Box
BF40401	135-9945	Mini Forceps with Spike	Alligator	Yellow	160	2	10/Box

## BIOPSY FORCEP FOR 2.8MM SCOPE CHANNELS STERILE ②

Order Number	Henry Schein Item Number	Description	Jaw Style	Sheath Color	Working Length (cm)	Minimum Channel Size (mm)	Package Units
BF40081	129-0179	Standard Forceps	Oval	Orange	180	2.8	10/Box
BF40091	129-0677	Standard Forceps	Oval	Blue	230	2.8	10/Box
BF40096	128-5668	Standard Forceps	Oval	Blue	230	2.8	50/Box
BF40111	129-0180	Standard Forceps with Spike	Oval	Orange	180	2.8	10/Box
BF40121	129-0678	Standard Forceps with Spike	Oval	Blue	230	2.8	10/Box
BF40126	128-5669	Standard Forceps with Spike	Oval	Blue	230	2.8	50/Box
BF40191	129-0184	Standard Forceps	Alligator	Orange	180	2.8	10/Box
BF40201	129-0679	Standard Forceps	Alligator	Blue	230	2.8	10/Box
BF40206	128-5670	Standard Forceps	Alligator	Blue	230	2.8	50/Box
BF40221	129-0185	Standard Forceps with Spike	Alligator	Orange	180	2.8	10/Box
BF40231	129-0680	Standard Forceps with Spike	Alligator	Blue	230	2.8	10/Box
BF40236	128-5671	Standard Forceps with Spike	Alligator	Blue	230	2.8	50/Box
BF40431	132-9020	Standard Forceps	Serrated	Orange	180	2.8	10/Box
BF40441	132-9017	Standard Forceps	Serrated	Blue	230	2.8	10/Box
BF40446	131-2212	Standard Forceps	Serrated	Blue	230	2.8	50/Box
BF40411	132-9018	Standard Forceps with Spike	Serrated	Orange	180	2.8	10/Box
BF40421	132-9019	Standard Forceps with Spike	Serrated	Blue	230	2.8	10/Box
BF40426	131-2211	Standard Forceps with Spike	Serrated	Blue	230	2.8	50/Box

## BIOPSY FORCEP FOR 3.2MM SCOPE CHANNELS STERILE ②

Order Number	Henry Schein Item Number	Description	Jaw Style	Sheath Color	Working Length (cm)	Minimum Channel Size (mm)	Package Units
BF40061	129-0177	Jumbo Biopsy	Oval	Blue	230	3.2	10/Box
BF40071	129-0178	Jumbo Biopsy with Spike	Oval	Blue	230	3.2	10/Box
BF40171	129-0182	Jumbo Biopsy	Alligator	Blue	230	3.2	10/Box
BF40181	129-0183	Jumbo Biopsy with Spike	Alligator	Blue	230	3.2	10/Box

## BIOPSY FORCEP HOT STERILE ②

Order Number	Henry Schein Item Number	Description	Jaw Style	Sheath Color	Working Length (cm)	Minimum Channel Size (mm)	Package Units
HF48521	129-0681	Biopsy Forceps Hot Lower	Oval	Red	230	2.8	10/Box
HF48541	131-2213	Biopsy Forceps Hot Lower	Serrated	Blue	230	2.8	10/Box

# PRECISION & CONTROL RETRIEVAL FORCEPS

Single-use grasping forceps for secure removal of tissue or foreign bodies.

- 3 Jaw types
- Ergonomic 3-ring handle
- Entry - Exit markers
- Sterile - Single Use Convenience



# DILATION

- ABC 3-STEP INFLATION BALLOONS
- MULTI-STAGE DILATION BALLOONS
- INFLATION DEVICE
- DILATION ACCESSORIES

## RETRIEVAL FORCEPS STERILE

Order Number	Henry Schein Item Number	Description	Jaw Opening (mm)	Sheath Color	Working Length (cm)	Minimum Channel Size (mm)	Package Units
GF61021	131-2207	Rat Tooth	6.3	Blue	230	2.8	10/Box
GF61041	131-2208	Alligator Tooth	8.1	Blue	230	2.8	10/Box
GF61061	131-2209	Mixed Tooth	8.1	Blue	230	2.8	10/Box






# ABC™ 3-STEP INFLATION BALLOONS

For a simple, affordable solution to balloon dilation, choose the ABC Dilation Balloon. Manufactured from a high-quality non-compliant nylon material, the ABC Balloon offers the radial strength required to accurately and consistently dilate strictures.

The balloon can be inflated in 3 steps, eliminating the need for multiple balloons in the same procedure. It offers an economical solution versus typical radial expansion balloons.

The balloon is available as a flexible fixed wire esophageal 8 cm balloon or a fixed wire esophageal, pyloric and colonic 5.5 cm balloon, both featuring rapid inflation/deflation performance.

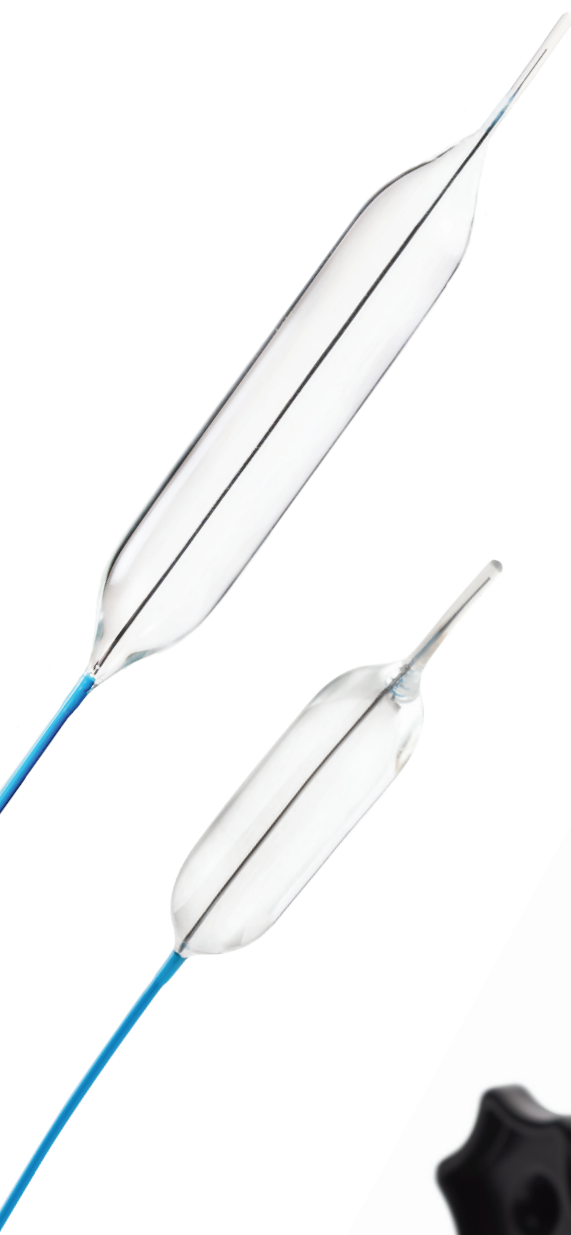

## KEY BENEFITS

### CHOICE OF MATERIAL

The durable nylon material delivers strength with accurate sizing. The non-compliant material provides the radial force required for consistent dilatation. Inflated in 3 steps, the ABC balloon provides an affordable alternative to radial expanding balloons.

### ROUNDED PROXIMAL SHOULDER

The combination of a rounded proximal shoulder and high-quality transparent nylon material facilitates visualization of the stricture during dilation.

### BALLOON LENGTH

Choose from a balloon length of 5.5 cm, which can be used as a universal device to dilate esophageal, pyloric or colonic strictures, or an 8 cm balloon depending on clinical need.

## ABC ESOPHAGEAL, PYLORIC, COLONIC DILATION BALLOON STERILE

Order Number	Henry Schein Item Number	Description	Balloon Diameter (mm)	Balloon Length (cm)	Working Length (cm)	Minimum Channel Size (mm)	Package Units
FW24165	139-6667	ABC Dilation Balloon	6	5.5	230	2.8	2/Box
FW24175	139-6668	ABC Dilation Balloon	8	5.5	230	2.8	2/Box
FW24185	139-6669	ABC Dilation Balloon	10	5.5	230	2.8	2/Box
FW24195	139-6670	ABC Dilation Balloon	12	5.5	230	2.8	2/Box
FW24205	139-6671	ABC Dilation Balloon	15	5.5	230	2.8	2/Box
FW24215	139-6672	ABC Dilation Balloon	18	5.5	230	2.8	2/Box
FW24225	139-6673	ABC Dilation Balloon	20	5.5	230	2.8	2/Box

## ABC ESOPHAGEAL DILATION BALLOON STERILE

Order Number	Henry Schein Item Number	Description	Balloon Diameter (mm)	Balloon Length (cm)	Working Length (cm)	Minimum Channel Size (mm)	Package Units
FW24075	128-5660	ABC Esophageal Dilation Balloon	6	8	180	2.8	2/Box
FW24085	128-5661	ABC Esophageal Dilation Balloon	8	8	180	2.8	2/Box
FW24095	128-5662	ABC Esophageal Dilation Balloon	10	8	180	2.8	2/Box
FW24105	128-5663	ABC Esophageal Dilation Balloon	12	8	180	2.8	2/Box
FW24125	128-5664	ABC Esophageal Dilation Balloon	15	8	180	2.8	2/Box
FW24145	128-5665	ABC Esophageal Dilation Balloon	18	8	180	2.8	2/Box
FW24155	128-5666	ABC Esophageal Dilation Balloon	20	8	180	2.8	2/Box

## INFLATION DEVICE STERILE

Order Number	Henry Schein Item Number	Description	Package Units
DB25011	131-4493	Micro-Tech Balloon Inflation Device 60 ml with lock	1/Box




# STRENGTH AND CONTROL MULTI-STAGE DILATION BALLOONS

Choose strength and control when dilating strictures with a dilation balloon.

Multi-Stage balloons feature a premium clear Pebax® balloon with rounded shoulders for balloon endoscopy

Fixed wire esophageal balloons offer fast, reliable performance

Pre-loaded .035 wire-guided balloons assist placement in tortuous or altered anatomy

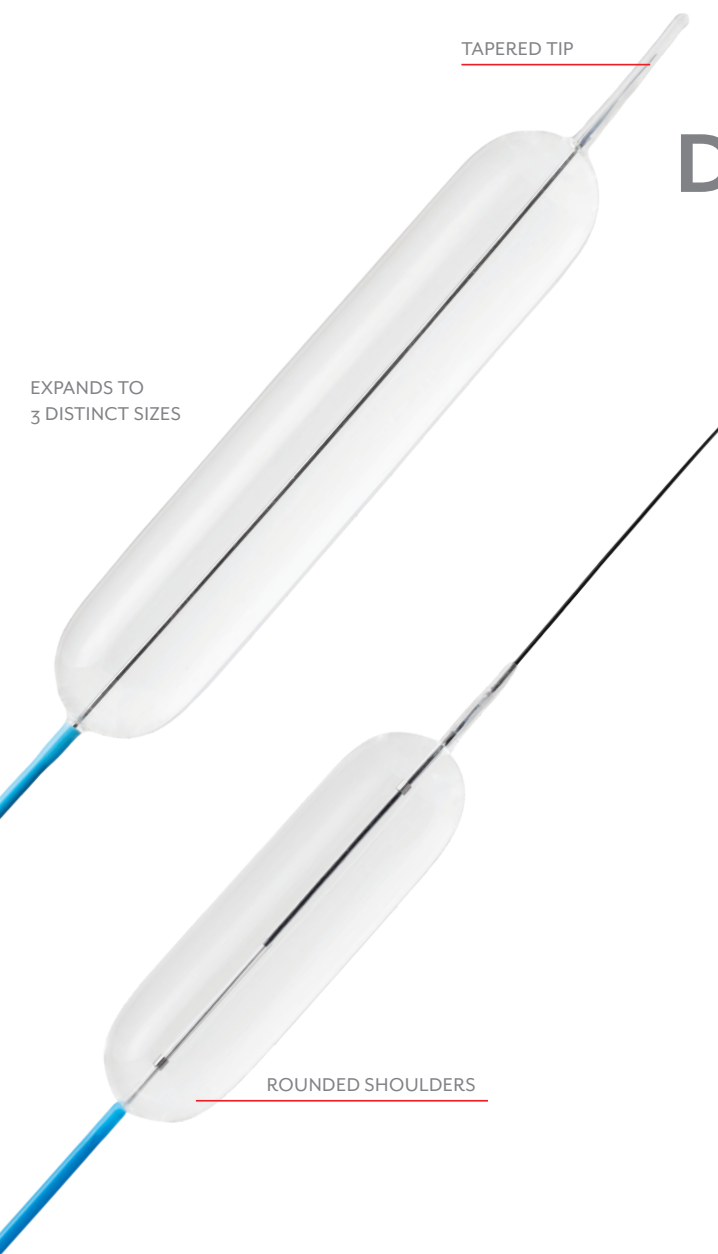

## WIRE-GUIDED ESOPHAGEAL, PYLORIC, COLONIC MULTI-STAGE BALLOONS STERILE

Order Number	Henry Schein Item Number	Description	Balloon Diameter (mm)	Balloon Length (cm)	Working Length (cm)	Minimum Channel Size (mm)	Package Units
MB23031	125-2734	Reliant™ Multi-Stage Dilatation Balloon	6-7-8	5.5	230	2.8	1/Box
MB23071	125-2736	Reliant™ Multi-Stage Dilatation Balloon	8-9-10	5.5	230	2.8	1/Box
MB23111	125-2738	Reliant™ Multi-Stage Dilatation Balloon	10-11-12	5.5	230	2.8	1/Box
MB23151	125-2740	Reliant™ Multi-Stage Dilatation Balloon	12-13.5-15	5.5	230	2.8	1/Box
MB23191	125-2742	Reliant™ Multi-Stage Dilatation Balloon	15-16.5-18	5.5	230	2.8	1/Box
MB23221	125-2744	Reliant™ Multi-Stage Dilatation Balloon	18-19-20	5.5	230	2.8	1/Box

## FIXED WIRE ESOPHAGEAL MULTI-STAGE BALLOONS STERILE

Order Number	Henry Schein Item Number	Description	Balloon Diameter (mm)	Balloon Length (cm)	Working Length (cm)	Minimum Channel Size (mm)	Package Units
FW24015	129-0163	Rapide™ Esophageal Dilatation Balloon	6-7-8	8	180	2.8	2/Box
FW24025	129-0164	Rapide™ Esophageal Dilatation Balloon	8-9-10	8	180	2.8	2/Box
FW24035	129-0165	Rapide™ Esophageal Dilatation Balloon	10-11-12	8	180	2.8	2/Box
FW24045	129-0166	Rapide™ Esophageal Dilatation Balloon	12-13.5-15	8	180	2.8	2/Box
FW24055	129-0167	Rapide™ Esophageal Dilatation Balloon	15-16.5-18	8	180	2.8	2/Box
FW24065	129-0168	Rapide™ Esophageal Dilatation Balloon	18-19-20	8	180	2.8	2/Box

## KEY BENEFITS

### CHOICE OF MATERIAL

The high quality Pebax material produces a balloon with the features and performance your providers demand. Accurate three-size radial expansion, rounded shoulders, clear visibility and strength with control.

### VISIBILITY UNDER FLUOROSCOPY

All wire-guided balloons feature two radiopaque markers at the balloon shoulder, allowing for accurate placement across the stricture. The dedicated inflation port permits the introduction of contrast to further enhance visibility if desired.

### TAPERED TIP

The rounded tip aids in the smooth passage of the balloon across the stricture. The wire-guided balloon is provided with a pre-loaded .035 guidewire that has a soft, flexible tip to assist in placement in tortuous or altered anatomy.



# SIMPLE SOLUTION BALLOON INFLATION DEVICE

The Micro-Tech Inflation Device was designed to accommodate smaller hands, featuring a lever-lock.



- Lever-lock controlled with the holding hand
- Plunger moved to prime, inflate, or create a vacuum
- Luminous dial includes ATMs on the outer scale for clear visualization
- Priced competitively
- 60 ml single-use syringe includes a fill mark for Micro-Tech Balloons for a significant value.



LOCKED



UNLOCKED

## INFLATION DEVICE STERILE


Order Number	Henry Schein Item Number	Description	Package Units
DB25011	131-4493	Micro-Tech Balloon Inflation Device 60 ml with lock	1/Box

# FREEDOM FROM REPROCESSING SINGLE USE GUIDEWIRE

Repeated use and reprocessing of reusable guidewires and cleaning brushes can result in worn, frayed, bent or otherwise damaged devices.

The use of disposable guidewires and disposable cleaning brushes can reduce the risk of device damage as well as the potential cross-contamination from using reprocessed accessories with reusable polyvinyl dilator systems.

The Single-Use Marked Spring Tip Guidewire and Dilator Cleaning Brush are for use with American and Savary dilators.

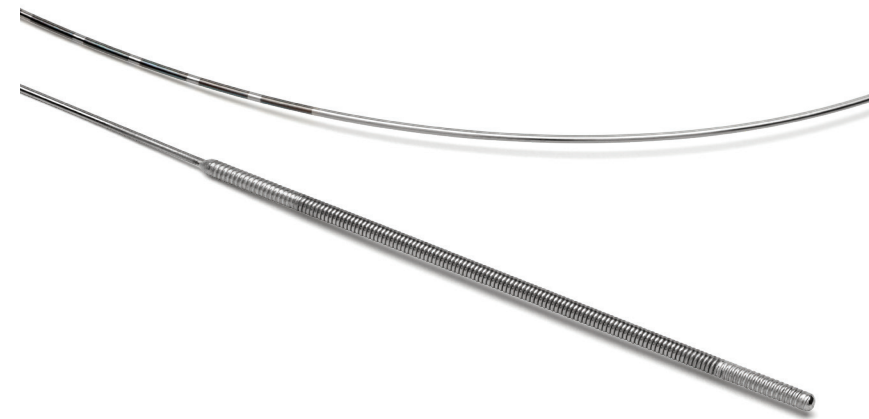

## MARKED SPRING TIP GUIDEWIRE STERILE

Order Number	Henry Schein Item Number	Description	Length (cm)	Package Units
DG25111	136-8998	Single-Use Marked Spring Tip Guide Wire for Esophageal Dilators	210	10/Box



## DILATOR CLEANING BRUSH

Order Number	Henry Schein Item Number	Description	Package Units
CB62411	135-6055	Dilator Cleaning Brush, 3mm for reusable American and Savary Dilators	10/Box

# HEMOSTASIS

---

INJECTION NEEDLES

LOCKADO HEMOSTASIS CLIP

SURECLIP HEMOSTASIS CLIP


HEMORRHOIDS ENDOSCOPIC  
MULTI-BAND LIGATORS

IN-SIGHT BANDER





# SHARP INNOVATION INJECTION NEEDLES

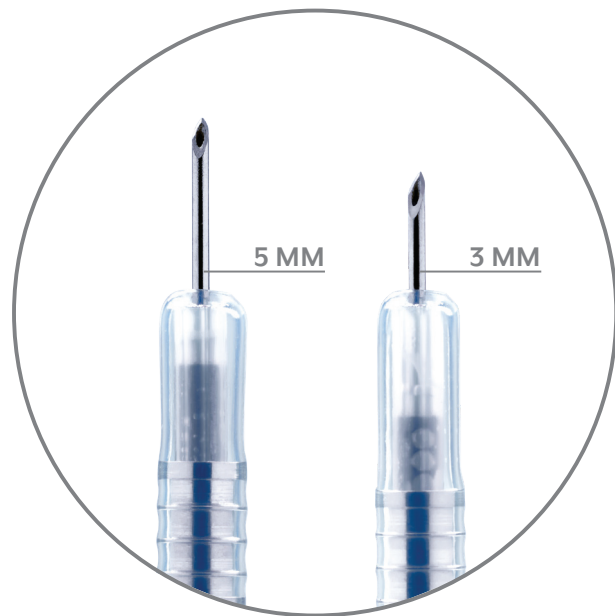

A simple, elegant solution, the 5 mm and 3 mm long lancet Injection Needles are offered in 22 and 25 gauge. Designed for single-handed operation, the innovative design limits accidental needle deployment during insertion and removal. The patented ratchet mechanism provides audible confirmation of deployment and minimizes needle bounce. A quick release button returns the needle fully into the sheath. The recessed metal end cap and firm catheter aid in tissue puncture and ensure needle performance through various scope positions.

## KEY BENEFITS

### PATENTED RATCHET DESIGN

Provides audible confirmation and secure locking of the deployed needle during use. As the needle is advanced, the ratchet gives tactile feedback and helps to reduce needle bounce.

### SHORT BEVEL TIP


Both the 5 mm and 3 mm needles feature a shorter three-bevel tip, designed to improve accuracy. The 3 mm PinPoint needle was designed to help target the submucosal space while the 5 mm allows for greater flexibility in therapeutic use.

### TARGET TISSUE DEPTH

With the shorter needle tip, PinPoint helps to consistently target the submucosal space. The 3 mm projection is sufficient to reach the submucosa in even the thickest intestinal tissue.

### ONE-HANDED OPERATION

Ergonomically designed one-handed operation enables confident needle deployment that enables rapid needle retraction.

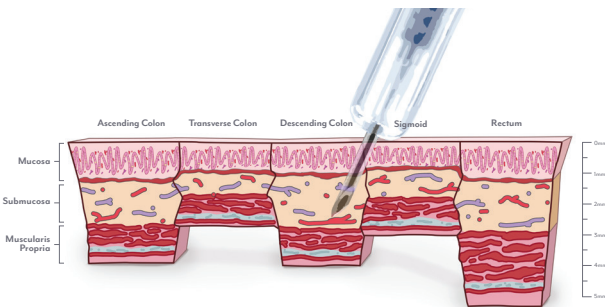

\*Tissue depth representative of thickest tissue recored in study, *Echoendoscopic characterization of the human colon*, Fernando M. Castro-Poças, Mário Dinis-Ribeiro, Tarcisio P. Araújo and Isabel Pedrotto

## INJECTION NEEDLES STERILE

Order Number	Henry Schein Item Number	Description	Needles Gauge/Length (mm)	Sheath Diameter (mm)	Working Length (cm)	Minimum Channel Size (mm)	Package Units
IN33131	129-0169	Injection Needle Upper	22/5	2.3	180	2.8	10/Box
IN33221	129-0171	Injection Needle Upper	25/5	2.3	180	2.8	10/Box
IN33151	129-0170	Injection Needle Lower	22/5	2.3	230	2.8	10/Box
IN33241	129-0172	Injection Needle Lower	25/5	2.3	230	2.8	10/Box
IN33301	139-5191	PinPoint Injection Needle Lower	22/3	2.3	230	2.8	10/Box
IN33331	139-5192	PinPoint Injection Needle Lower	25/3	2.3	230	2.8	10/Box

# LOCK ONTO TISSUE LOCKADO™ HEMOSTASIS CLIP



Hemostasis clips need to provide accuracy, reliability and precise functionality. Lockado allows unrestricted opening and closing so it can be repositioned until the desired position is achieved.

With two sets of opposing atraumatic teeth in the arms of the clip, Lockado delivers a reliable and secure grip for consistent performance. Driven by Micro-Tech Endoscopy's proven clip delivery system, the Lockado clip provides exceptional reliability and tissue management.

- Jaw style works for multiple clinical applications
- Can be used with side-viewing endoscope
- Shorter stem make clips less obtrusive
- Variety of sizes: 11 mm, 16 mm and NEW 22 mm arm opening widths
- Indicated for anchoring jejunal feeding tubes

## KEY BENEFITS

### LOCKING SPOOL MECHANISM

Mechanism provides tactile and audible feedback of clip deployment. Spool is locked in place and clip is released automatically.

### RIBBED CATHETER

The lubricious ribbed catheter surrounds a torque wire that provides 1:1 response and 360° rotation without inadvertent clip deployment or helicoptering.

### ROTATION GRIP

Secure grip allows for 360° bidirectional controlled rotation and placement of clip on targeted tissue.

## LOCKADO HEMOSTASIS CLIP, REPOSITIONABLE, ROTATABLE STERILE

Order Number	Henry Schein Item Number	Opening Width (mm)	Sheath Diameter (mm)	Working Length (cm)	Minimum Channel Size (mm)	Packing Units
LO30861	138-6525	11	2.6	235	2.8	10/Box
LO30865	138-6526	11	2.6	235	2.8	2/Box
LO30891	138-6527	16	2.6	235	2.8	10/Box
LO30895	138-6529	16	2.6	235	2.8	2/Box
LO35041	TBD	22	2.6	235	2.8	10/Box
LO35045	TBD	22	2.6	235	2.8	2/Box




# CLIP WITH CONFIDENCE SURECLIP®

Clips need to be reliable. They need to be accurate. They need to allow you the flexibility to reposition or rotate as much as is required to deliver better outcomes.


Accurate positioning prior to deployment can reduce both procedure time and the number of clips needed to achieve hemostasis. SureClip achieves this by design delivering outstanding repositionability and reliable rotation prior to deployment, in various scope positions. SureClip's short stem aids placement in narrow lumen and improves visibility.

The proprietary clip design provides reliable deployment, may improve retention, and offers a choice of jaw sizes.

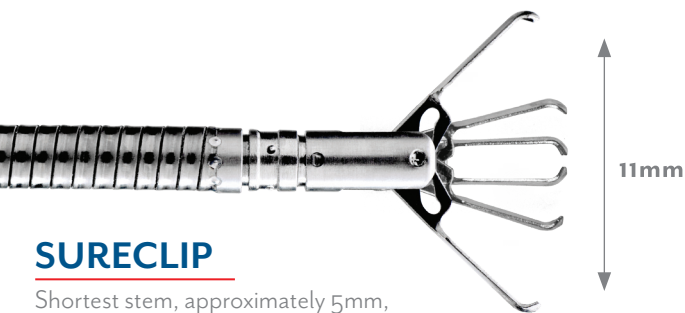

## SURECLIP

Shortest stem, approximately 5mm, for use in narrow lumen

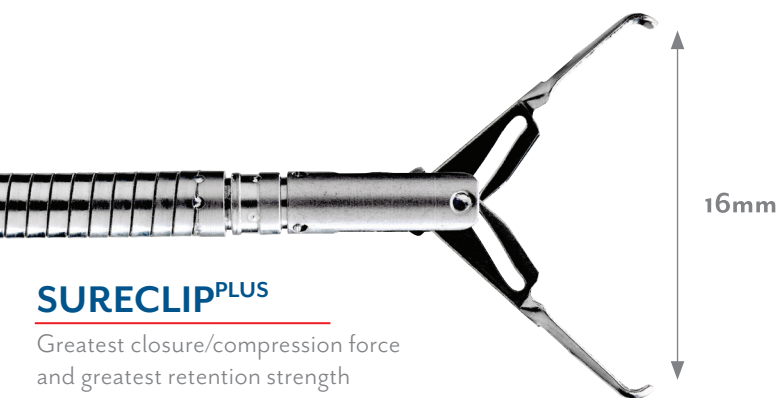


## SURECLIPPLUS

Greatest closure/compression force and greatest retention strength



## SURECLIPMAX

Widest jaw opening to accommodate greater tissue volume

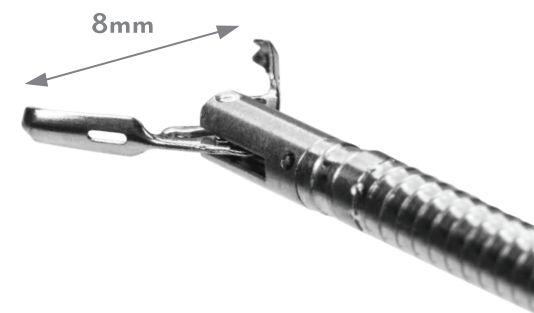

## SURECLIPMINI

Lowest jaw profile for unobstructed view during placement

## KEY BENEFITS

### REPOSITIONING

SureClip's unique design permits opening and closing the jaw prior to deployment. Being able to reposition a clip may help improve placement accuracy. Fenestrations on the SureClipPlus, SureClipmini and SureClipmax accommodate tissue and may enhance retention.

### RELIABLE ROTATION

SureClip can be rotated, helping to provide the correct orientation for tissue approximation or defect closure. The rotation handle on the SureClip improves performance and enhances the user experience.

### SHORT STEM

A shorter stem makes the clip less obtrusive, improving visualization of the target area, particularly when multiple clips are placed in close proximity.



### SURECLIP<sup>MINI</sup> LOWEST PROFILE STERILE

Order Number	Henry Schein Item Number	Opening Width (mm)	Sheath Diameter (mm)	Working Length (cm)	Minimum Channel Size(mm)	Package Units
RC30415	132-5180	8	Max 2.6	235	2.8	2/Box
RC30411	132-5187	8	Max 2.6	235	2.8	10/Box



### SURECLIP SHORTEST STEM STERILE

Order Number	Henry Schein Item Number	Opening Width (mm)	Sheath Diameter (mm)	Working Length (cm)	Minimum Channel Size(mm)	Package Units
RC30445	132-5724	11	Max 2.6	235	2.8	2/Box
RC30441	132-5723	11	Max 2.6	235	2.8	10/Box



### SURECLIP<sup>PLUS</sup> GREATEST CLOSURE/COMPRESSION FORCE STERILE

Order Number	Henry Schein Item Number	Opening Width (mm)	Sheath Diameter (mm)	Working Length (cm)	Minimum Channel Size(mm)	Package Units
RC30385	128-5657	16	Max 2.6	235	2.8	2/Box
RC30381	128-5655	16	Max 2.6	235	2.8	10/Box



### SURECLIP<sup>MAX</sup> WIDEST JAW OPENING\* STERILE

Order Number	Henry Schein Item Number	Opening Width (mm)	Sheath Diameter (mm)	Working Length (cm)	Minimum Channel Size(mm)	Package Units
RC30545	1366309	17	Max 2.6	235	2.8	2/Box
RC30541	1366310	17	Max 2.6	235	2.8	10/Box

\*Jaws may not open to maximum reach until placed in contact with tissue





## EFFECTIVE LIGATION MULTIPLE BAND LIGATORS



Endoscopic banding has been demonstrated as an effective and simple treatment for internal hemorrhoids and varices. Hemorrhoid banding is readily performed in an ambulatory setting and can result in less pain and have a shorter recovery period than a surgical hemorrhoidectomy.

The Multi-Band Ligator Sets from Micro-Tech offer cost-effective solutions for the treatment of hemorrhoids and varices during routine endoscopy procedures.

### KEY BENEFITS

#### EXCELLENT VISIBILITY

A clear ligating barrel with rear mounted bands optimizes the endoscopic view while maintaining suction.

#### BAND DESIGN

The square band design provides consistent shape and reliable retention. The hemorrhoid set comes with four pre-loaded latex-free bands for routine colonoscopies. The upper ligator comes pre-loaded with six latex-free bands and can be used in upper and lower endoscopic procedures.

#### USER-FRIENDLY SYSTEM

The handle unit tie hole enables quick and easy loading of the pull wire. The alert band notifies the operator when only one band remains.

#### MULTI-BAND LIGATOR STERILE

Order Number	Henry Schein Item Number	Description	Minimum Working Channel (mm)	Working Length (cm)	Package Units
HL32161	131-9042	Endoscopic Hemorrhoid Multi-Band Ligator Set, 4 Latex Free Bands, compatible with 9.4-13 mm scope diameter	2.8	220	1/Box
ML32651	138-5348	Endoscopic Multi-Band Ligator Set, 6 Latex Free Bands, compatible with 9.4-13 mm scope diameter	2.8	145	1/Box

## IN-OFFICE LIGATION IN-SIGHT™ BANDER



Hemorrhoid banding has been demonstrated as an effective and simple treatment for internal hemorrhoids, and can be readily performed in an office or ambulatory setting. Banding can also result in less pain with a shorter recovery period than a surgical hemorrhoidectomy.

The In-Sight Multi-Band Ligator is the only hand-held ligation device with its own light source – providing illumination, suction, direct visualization and ligation capabilities in a single-use device.

- Pre-loaded with 4 ligator bands
- Comfortable and ergonomic design for easy ligation and patient comfort
- Unique LED system maximizes visualization safely and efficiently



#### IN-OFFICE MULTI-BAND LIGATOR

Order Number	Henry Schein Item Number	Description	Package Units
QL30011	138-8285	In-Sight Multi-Band Ligator	10/Box

# INFECTION PREVENTION

JACKET DISPOSABLE BIOPSY VALVE

CATCHALL TRANSPORT PAD

CLEANING BRUSHES

ACCESSORY KITS

## JACKET™ DISPOSABLE BIOPSY VALVE

The Jacket Disposable Biopsy Valve is used to cover the opening to the biopsy/suction channel of gastrointestinal endoscopes.

Its unique design offers improved valve performance by minimizing leakage and helping to preserve insufflation to maintain adequate visualization.

- Double valve system to minimize leakage and prevent spraying
- Improved seal with three points of contact
- Internal funnel seal hugs and flexes with inserted accessories without gaping
- Compatible with Olympus, Pentax and Fuji GI endoscopes



### JACKET DISPOSABLE BIOPSY VALVE ②

Order Number	Henry Schein Item Number	Description	Package Units
QL70021	138-6532	Jacket Disposable Biopsy Valve, Olympus/Fujinon, Non-Sterile (Blue)	100/Box
QL70051	138-6534	Jacket Disposable Biopsy Valve, Pentax Medical, Non-Sterile (Clear)	100/Box



# PROTECT, COLLECT, CLEAN GI ACCESSORIES



## CATCHALL SCOPE TRANSPORT PAD <sup>Ⓢ</sup>

The CatchAll Transport Pad offers compliance and protection while transporting contaminated scopes.

Order Number	Henry Schein Item Number	Description	Package Units
CP62121	131-2218	CatchAll Scope Transport Pad	25/Box



## CLEANING BRUSHES <sup>Ⓢ</sup>

Micro-Tech's Cleaning Brushes combine performance and compliance. Available in combination kits, the brush range features 6/10mm nylon brush heads designed for effective gentle cleaning.

Order Number	Henry Schein Item Number	Description	Package Units
QL22011	131-2220	Single ended 6mm channel brush	50/Box
QL22021	131-2222	Double ended 6mm channel brush	50/Box
QL22031	131-2221	Single ended 6mm channel and attached 10mm valve brush	50/Box
QL22041	131-2215	Double ended 5mm-10mm valve brush	50/Box
QL22051	131-2216	Single ended 6mm channel brush + 5mm-10mm valve brush	50/Box
QL22061	131-2217	Double ended 6mm channel brush + 5mm-10mm valve brush	50/Box

# CONVENIENT COMPLIANCE GI ACCESSORY KITS



Procedure kits provide the essential supplies needed for most cases in one convenient package. Using preconfigured accessory kits helps to standardize processes, reduce product handling, and simplify supply ordering.

Choose from four different kits that provide the items required for nearly every endoscopic procedure. By bundling routine items into standardized kits, the number of touchpoints, and the potential for cross-contamination may be reduced.

- Basic Kit
- Upper Kit
- Set-Up Kit
- Compliance Kit

	Catch-All Transport Pad	Biohazard Label	Bowl	Sponge (w/ detergent)	Gauze (10)	Lubricating Gel (1.4 oz)	Cleaning Brush	Bite Block	Suction Tubing
Basic Kit	•	•	•	•	•	•			
Upper Kit	•	•	•	•	•	•		•	
Set-Up Kit	•	•	•	•	•	•	•		•
Compliance Kit	•	•	•	•					

## GI ACCESSORY KITS <sup>Ⓢ</sup>

Order Number	Henry Schein Item Number	Description	Package Units
QL22531	140-0557	Set-Up Endoscopy Kit	12/Box
QL22541	140-0704	Compliance Endoscopy Kit	12/Box
QL22511	140-0706	Basic Endoscopy Kit	12/Box
QL22521	140-0705	Upper Endoscopy Kit	12/Box
QL22711	140-1637	Lubricating Gel, 1.4 oz Pouch	100/Box
QL22811	140-1682	Enzymatic Sponge	50/Box





# POLYPECTOMY

ENDOINK TATTOO

HOT SNARES

COLD SNARES

BLUE BEACON GI MARKER

ENTRAP POLYP RETRIEVAL NET

## ENDOSCOPIC TATTOO ENDOINK™

EndoINK is a single-use fluid designed to endoscopically mark lesions in the gastrointestinal tract. Developed with a suspension liquid containing highly purified carbon particles, EndoINK creates a marker to facilitate surgical removal and examination.



### KEY BENEFITS

**RELIABLE MARKER**

Endoscopic tattoo for reliable marking of lesions and areas of reference in the GI tract

**C 10ML SYRINGE**

Conveniently available in 10 ml configuration to accommodate catheter priming

**EASY TO USE**

Sterile fluid supplied in a capped 10 ml syringe for single patient use

**ENDOINK** STERILE

Order Number	Henry Schein Item Number	Description	Package Units
NK62721	139-2221	EndoINK Endoscopic Marker, 10ml, Sterile Fluid	5/Box




# PRECISE CUT HOT SNARES

Our line of hot and cold snares features effective cutting wires and user-friendly handles at cost-effective prices.

With the innovative power connector, the snare can be operated and rotated with ease

Our hot snares are available in a compliant braided wire or a firmer wire for effective tissue acquisition and resection

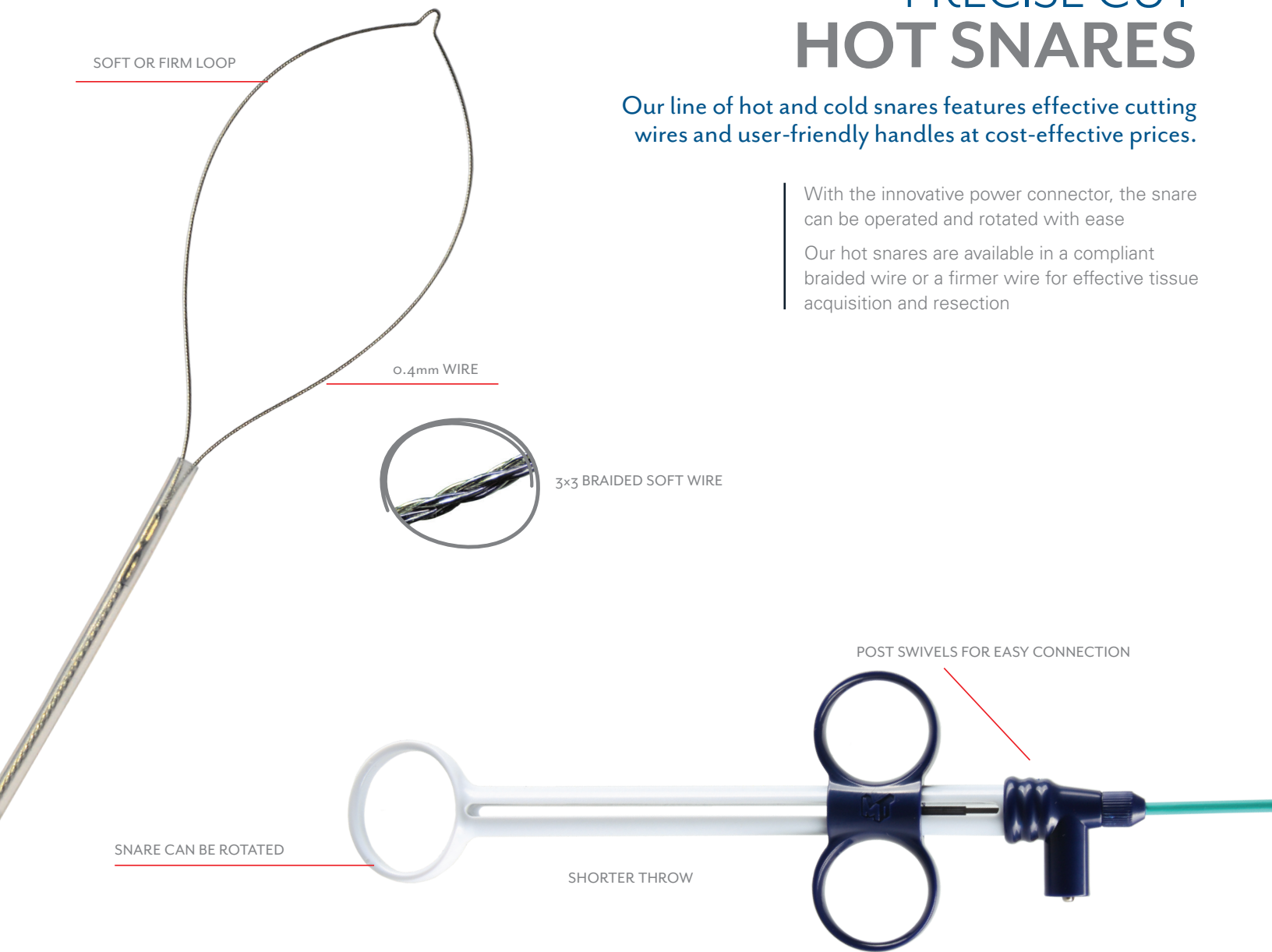

## HOT SNARES BRAIDED STERILE

Order Number	Henry Schein Item Number	Description	Loop Diameter (mm)	Sheath Diameter (Fr)	Working Length (cm)	Minimum Channel Size (mm)	Package Units
PS51011	129-0186	Mini Braided Oval	10	7	230	2.8	10/Box
PS51021	129-0673	Small Braided Oval	15	7	230	2.8	10/Box
PS51031	129-0674	Medium Braided Oval	20	7	230	2.8	10/Box
PS51041	129-0675	Large Braided Oval	24	7	230	2.8	10/Box
PS51051	129-0187	Extra Large Braided Oval	30	7	230	2.8	10/Box

## HOT SNARES FIRM STERILE

Order Number	Henry Schein Item Number	Description	Loop Diameter (mm)	Sheath Diameter (Fr)	Working Length (cm)	Minimum Channel Size (mm)	Package Units
PS51071	129-0188	Mini Oval Firm	10	7	230	2.8	10/Box
PS51081	129-0189	Small Oval Firm	15	7	230	2.8	10/Box
PS51091	129-0190	Medium Oval Firm	20	7	230	2.8	10/Box
PS51101	129-0191	Large Oval Firm	24	7	230	2.8	10/Box
PS51111	129-0192	Extra Large Oval Firm	30	7	230	2.8	10/Box
PS51121	129-0193	Jumbo Oval Firm	36	7	230	2.8	10/Box

## HOT SNARES HEX FIRM STERILE

Order Number	Henry Schein Item Number	Description	Loop Diameter (mm)	Sheath Diameter (Fr)	Working Length (cm)	Minimum Channel Size (mm)	Package Units
PS52031	129-0194	Small Hex Firm	15	7	230	2.8	10/Box
PS52041	129-0195	Large Hex Firm	25	7	230	2.8	10/Box

## KEY BENEFITS

### SHORTER THROW HANDLE

The handle throw is adjusted according to the snare size. Small and Medium snares have a shorter throw than Large or Extra Large snares. A shorter throw improves comfort for users with smaller hands.

### BRAIDED WIRE

The braided wire combines the tissue acquisition features of a firm snare with the flexibility and smooth deployment of a soft snare. Suitable for daily use in facilities looking to minimize inventory.

### INNOVATIVE POWER CONNECTION


The handle features a dedicated power connection, separate from the finger ring slider. This feature permits unrestricted movement and rotation of the snare handle and may reduce procedural delays.




OVAL OR DIAMOND LOOP

## EFFICIENT RESECTION COLD SNARES

- Designed specifically to resect diminutive polyps
- Available in two sizes – 10 mm and 15 mm
- Thin cutting wire provides for a guillotine-type cut
- Loop shape assists in capturing smaller polyps
- Stiffer catheter aids in insertion and provides support during resection

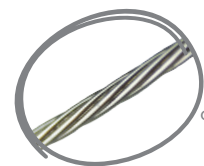

0.23mm WIRE

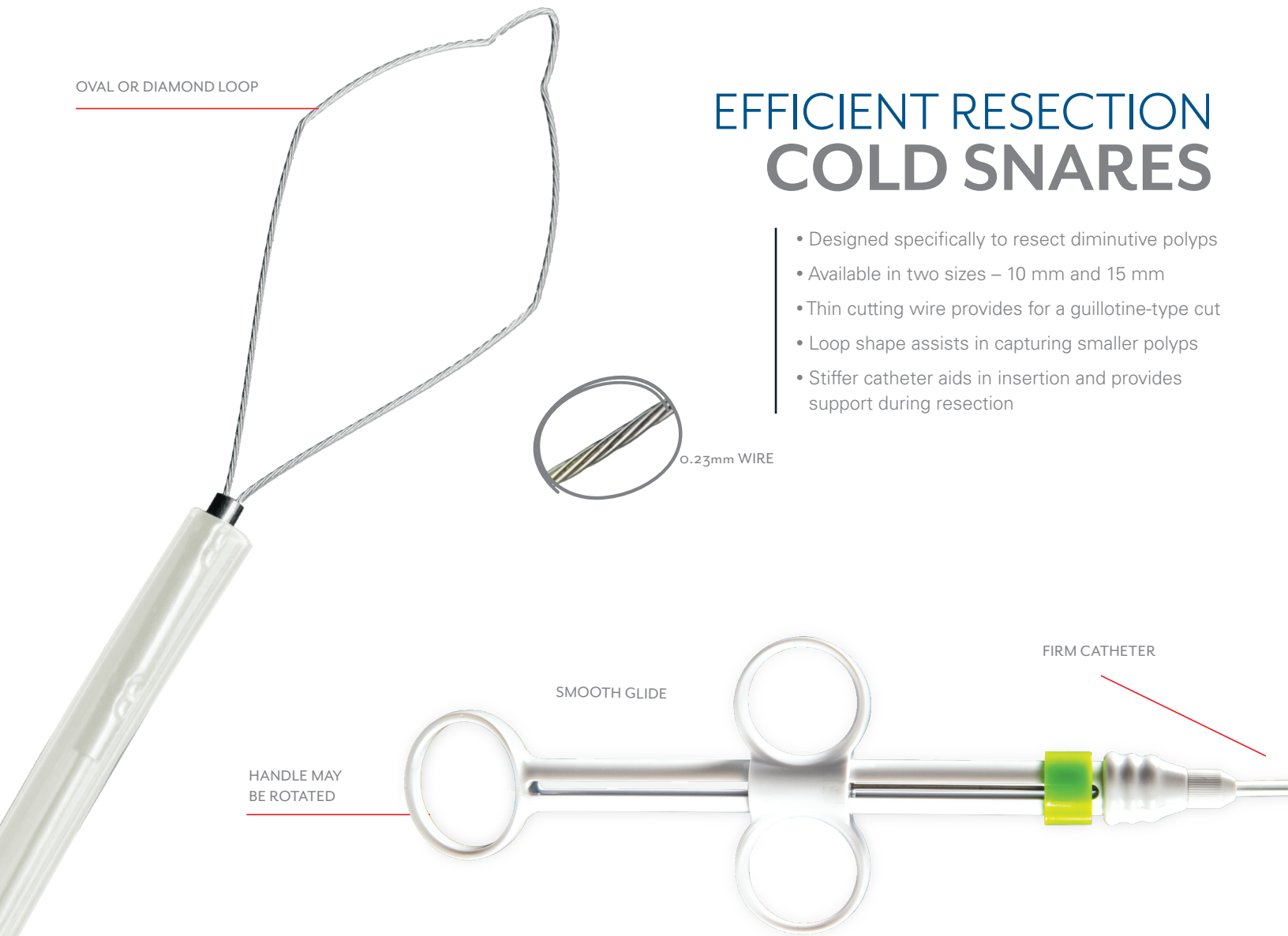

### KEY BENEFITS

#### DEDICATED COLD SNARE

Specifically designed for resecting diminutive polyps using a cold snare technique, the oval or diamond shape loop and flexible wire, capture smaller polyps with ease.

#### THIN WIRE

The thin 0.23 mm cutting wire provides a guillotine cut that reduces the force needed to complete the resection and may limit tearing or snagging.

#### COMFORTABLE


Accommodating various hand sizes, the ergonomic three-ring handle features a shortened throw and comfortable grip. The stiff catheter provides support for a smooth effortless cut.

### COLD SNARES STERILE

Order Number	Henry Schein Item Number	Description	Loop Diameter (mm)	Sheath Diameter (Fr)	Working Length (cm)	Minimum Channel Size (mm)	Package Units
CS50011	129-0671	Diamond Thin Wire	10	7	230	2.8	10/Box
CS50021	129-0672	Diamond Thin Wire	15	7	230	2.8	10/Box
CS50031	129-0160	Oval Thin Wire	10	7	230	2.8	10/Box
CS50041	129-0161	Oval Thin Wire	15	7	230	2.8	10/Box

## TRUE BLUE BLUE BEACON™ G.I. MARKER

- Useful in GI procedures that require a contrast stain to temporarily mark or highlight tissue
- Can be diluted and introduced to the irrigation, sprayed or injected around the area of interest
- The 0.4% Indigo Carmine solution is supplied sterile in a convenient luer lock 5ml syringe, 5 syringes per box.
- The 0.8% Indigo Carmine solution is supplied sterile in a 10ml single use vial, 10 vials per box.

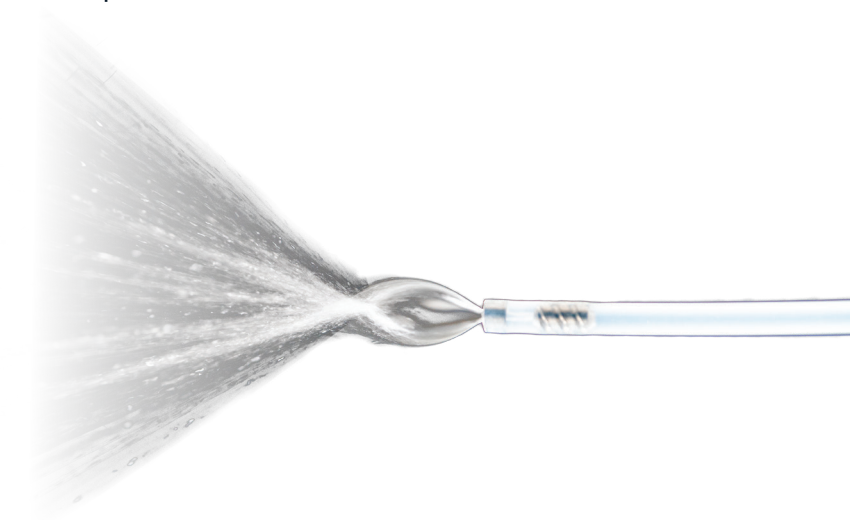

### BLUE BEACON G.I. MARKER STERILE

Order Number	Henry Schein Item Number	Description	Package Units
IC62021	128-5667	Blue Beacon GI Marker 5ml Syringe Indigo Carmine 0.4% Sterile	5/Box
IC62061	137-5185	Blue Beacon GI Marker 10ml Single Use Vial Indigo Carmine 0.8% Sterile	10/Box

This solution is NOT FOR INTRAVENOUS OR INTRAMUSCULAR INJECTION. It is a single use device indicated for temporary marking of the GI tract.


### SPRAY CATHETER STERILE

Order Number	Henry Schein Item Number	Description	Working Length (cm)	Minimum Channel Size (mm)	Package Units
SC62241	133-5544	Spray Catheter	220	2.2	20/Box
SC62301	133-5545	Spray Catheter	250	2.8	20/Box



# DURABLE AND GENTLE ENTRAP™ POLYP RETRIEVAL NET

The Entrap retrieval system delivers a quick and efficient retrieval solution for polyps and tissue for histopathologic review. Utilizing advanced materials and production techniques, the Entrap facilitates improved consistency and control of the net deployment, position, and closure.

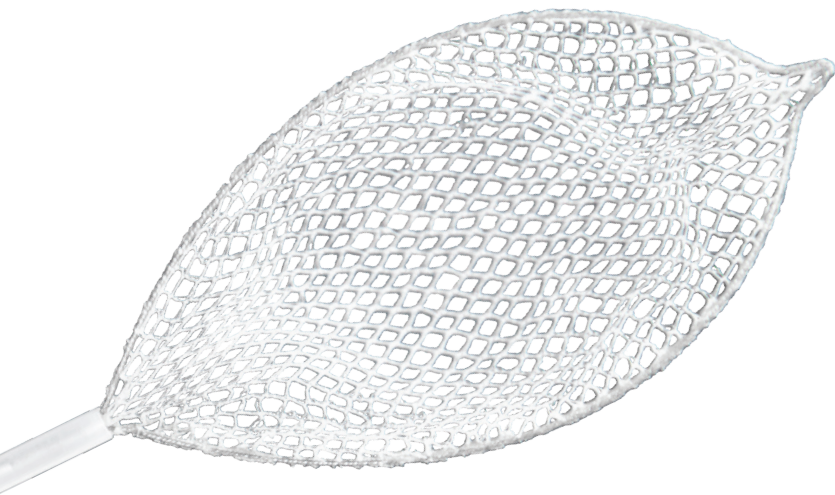

## KEY BENEFITS

### REFINED NET

Improved visibility through the net material for identification of the retrieval target. Durable enough for repeated deployment, but soft to gently capture and preserve the specimen.

### ADVANCED PROCESSES

Cinching, tearing, and binding of the net material is reduced by the minimal net selvage and bonded net attachment. This increases efficiency and minimizes the time required to reliably capture specimens.

### SMOOTH AND STRONG

A stiffer catheter material improves device insertion into the endoscope and positioning of the net. The additional columnar strength allows for gentle movement of a specimen for optimal capture.

## ENTRAP RETRIEVAL SYSTEM

Order Number	Henry Schein Item Number	Description	Loop Diameter (mm)	Sheath Diameter (Fr)	Working Length (cm)	Minimum Channel Size (mm)	Package Units
NT54511	133-5543	Entrap Polyp Retrieval Net	30	7	230	2.8	10/Box
NT54515	135-6057	Entrap Polyp Retrieval Net	30	7	230	2.8	2/Box

# GI ACCESSORIES

BREATHE EZ BITE BLOCK

BITE BLOCKS

POLYP TRAPS

PROTRAP LUMA



# PREMIUM MATERIAL AND ENHANCED DESIGN BREATHE EZ BITE BLOCK

Micro-Tech Bite Blocks offer comfort and protection during upper endoscopic procedures. The new Breathe EZ Bite Block was designed with a special, softer bite plate that helps to keep the lower jaw forward and the airway open.

- Unique design that may assist in improving oxygenation during endoscopic procedures
- Rigid strap to prevent the device from moving out of position
- Softer overlay that helps protect patient's teeth



## BREATHE-EZ BITE BLOCK

Order Number	Henry Schein Item Number	Description	Package Units
QL26021	TBD	Breathe-EZ Bite Block, Adult, 60Fr, with strap, Single Use	50/Box

# GI ACCESSORIES

Micro-Tech Bite Blocks offer comfort and protection during upper endoscopic procedures.



## BITE BLOCKS

Micro-Tech Bite Blocks offer comfort and protection during upper endoscopic procedures. Available in 38Fr, 60Fr and 60Plus sizes, with or without dental rim.

Order Number	Henry Schein Item Number	Description	Package Units
BB60011	125-2687	Bite Block, Adult, With Strap	100/Box
BB60021	125-2688	Bite Block, Pediatric, With Strap	100/Box
BB60041	131-2210	Bite Block, Adult w/o dental rim. 60+Fr with strap	100/Box
BB60051	133-5546	Bite Block, Adult, 60+ Fr with strap	100/Box



## POLYP TRAPS

Polyp Retrieval Traps offer improved efficiency, reduced handling and cost savings. The unique Quad trap features removable collection trays and a secondary filter for peace of mind.

Order Number	Henry Schein Item Number	Description	Package Units
QL54012	136-8983	In-Line Suction Polyp Trap	50/Box
PC54011	129-0676	Polyp Trap With 4 Chambers	25/Box
PC54021	131-2219	Polyp Trap single chamber	25/Box
PC54031	131-2214	Quad Polyp Trap with removable chambers	25/Box



# CLEARLY, THE ULTIMATE POLYP TRAP PROTRAP LUMA



The ProTrap Luma™, in-line polyp trap and specimen tray consist of a collection hose, a trap body, a tray, and a glow stick. Pro-Trap's unique design offers the potential for touch-less sample retrieval. The specimen tray can be removed and placed directly into the specimen jar, minimizing handling and the associated risks. Additional trays are available as needed, and the trap is complemented with a light source for ease of use.

- Bright and long-lasting glow stick helps aid in retrieval visualization
- The entire tray fits inside the specimen jar
- Built-in magnifying window for easy visualization of specimen



## PROTRAP LUMA ⓘ

Order Number	Henry Schein Item Number	Description	Package Units
TP54411	TBD	ProTrap Luma In-Line Polyp Trap	25/Box
TP54421	TBD	ProTrap Luma In-Line Polyp Trap Specimen Tray	50/Box

# PRODUCT REFERENCE CHART

Order Number	Ref #	UPN Pouch	Henry Schein Item Number	Product Family
BB60011	AC01-103.A	BB60010	125-2687	GI Accessories
BB60021	AC01-103.P	BB60020	125-2688	GI Accessories
BB60041	AC01-108.M	BB60040	131-2210	GI Accessories
BB60051	AC01-109.M	BB60050	133-5546	GI Accessories
BF40021	NBF02-11018160	BF40020	129-0175	Biopsy
BF40051	NBF02-11118160	BF40050	129-0176	Biopsy
BF40061	NBF03-01030230	BF40060	129-0177	Biopsy
BF40071	NBF03-01130230	BF40070	129-0178	Biopsy
BF40081	NBF03-11023180	BF40080	129-0179	Biopsy
BF40091	NBF03-11023230	BF40090	129-0677	Biopsy
BF40096	NBF03-11023230	BF40094	128-5668	Biopsy
BF40111	NBF03-11123180	BF40110	129-0180	Biopsy
BF40121	NBF03-11123230	BF40120	129-0678	Biopsy
BF40126	NBF03-11123230	BF40124	128-5669	Biopsy
BF40141	NBF12-11018160	BF40140	129-0181	Biopsy
BF40171	NBF13-01030230	BF40170	129-0182	Biopsy
BF40181	NBF13-01130230	BF40180	129-0183	Biopsy
BF40191	NBF13-11023180	BF40190	129-0184	Biopsy
BF40201	NBF13-11023230	BF40200	129-0679	Biopsy
BF40206	NBF13-11023230	BF40204	128-5670	Biopsy
BF40221	NBF13-11123180	BF40220	129-0185	Biopsy
BF40231	NBF13-11123230	BF40230	129-0680	Biopsy
BF40236	NBF13-11123230	BF40234	128-5671	Biopsy
BF40286	LCNBF13-10024230	BF40284	133-4667	Biopsy
BF40306	LCNBF13-10124230	BF40304	133-4666	Biopsy
BF40401	NBF12-11118160	BF40400	135-9945	Biopsy
BF40411	NBF63-11123180	BF40410	132-9018	Biopsy
BF40421	NBF63-11123230	BF40420	132-9019	Biopsy
BF40426	NBF63-11123230	BF40424	131-2211	Biopsy
BF40431	NBF63-11023180	BF40430	132-9020	Biopsy
BF40441	NBF63-11023230	BF40440	132-9017	Biopsy
BF40446	NBF63-11023230	BF40444	131-2212	Biopsy
CB62411	CB10-0030/09-S	CB62410	135-6055	Dilation
CP62121	CTP-01025	CP62120	131-2218	Infection Prevention
CS50011	CS3-11023230	CS50010	129-0671	Polypectomy
CS50021	CS3-11523230	CS50020	129-0672	Polypectomy
CS50031	CS3-21023230	CS50030	129-0160	Polypectomy
CS50041	CS3-21523230	CS50040	129-0161	Polypectomy
CY49021	CYB-20S-18C	CY49020	133-5541	Biopsy
CY49051	CYB-30S-24C	CY49050	133-5542	Biopsy
DB25011	BID-15-60-01	DB25010	131-4493	Dilation
DB25011	BID-15-60-01	DB25010	131-4493	Dilation
DG25111	DG-18-08-21	DG25110	136-8998	Dilation
FW24015	FWB-06/80-6/18-M	FW24014	129-0163	Dilation
FW24025	FWB-08/80-6/18-M	FW24024	129-0164	Dilation
FW24035	FWB-10/80-6/18-M	FW24034	129-0165	Dilation





Order Number	Ref #	UPN Pouch	Henry Schein Item Number	Product Family
FW24045	FWB-12/80-6/18-M	FW24044	129-0166	Dilation
FW24055	FWB-15/80-6/18-M	FW24054	129-0167	Dilation
FW24065	FWB-18/80-6/18-M	FW24064	129-0168	Dilation
FW24075	FWB-06/80-6/18-S	FW24074	128-5660	Dilation
FW24085	FWB-08/80-6/18-S	FW24084	128-5661	Dilation
FW24095	FWB-10/80-6/18-S	FW24094	128-5662	Dilation
FW24105	FWB-12/80-6/18-S	FW24104	128-5663	Dilation
FW24125	FWB-15/80-6/18-S	FW24124	128-5664	Dilation
FW24145	FWB-18/80-6/18-S	FW24144	128-5665	Dilation
FW24155	FWB-20/80-6/18-S	FW24154	128-5666	Dilation
FW24165	FWB-06/55-6/23-S	FW24164	139-6667	Dilation
FW24175	FWB-08/55-6/23-S	FW24174	139-6668	Dilation
FW24185	FWB-10/55-6/23-S	FW24184	139-6669	Dilation
FW24195	FWB-12/55-6/23-S	FW24194	139-6670	Dilation
FW24205	FWB-15/55-6/23-S	FW24204	139-6671	Dilation
FW24215	FWB-18/55-6/23-S	FW24214	139-6672	Dilation
FW24225	FWB-20/55-6/23-S	FW24224	139-6673	Dilation
GF61021	NGF03-11023230	GF61020	131-2207	Biopsy
GF61041	NGF13-11023230	GF61040	131-2208	Biopsy
GF61061	NGF33-11023230	GF61060	131-2209	Biopsy
HF48521	HBF55-11023230	HF48520	129-0681	Biopsy
HF48541	HBF65-11023230	HF48540	131-2213	Biopsy
HL32161	HMBLS-4F-NL	HL32160	131-9042	Hemostasis
IC62021	P-DYZ-04-05	IC62020	128-5667	Polypectomy
IC62061	DYZ-08-10	IC62060	137-5185	Polypectomy
IN33131	IN12-225231802	IN33130	129-0169	Hemostasis
IN33151	IN12-225232302	IN33150	129-0170	Hemostasis
IN33221	IN12-255231802	IN33220	129-0171	Hemostasis
IN33241	IN12-255232302	IN33240	129-0172	Hemostasis
IN33301	IN12-223232302	IN33300	139-5191	Hemostasis
IN33331	IN12-253232302	IN33330	139-5192	Hemostasis
LO30861	LOCK-D-26-235-C-N	LO30860	138-6525	Hemostasis
LO30865	LOCK-D-26-235-C-N	LO30860	138-6526	Hemostasis
LO30891	LOCK-F-26-235-C-N	LO30890	138-6527	Hemostasis
LO30895	LOCK-F-26-235-C-N	LO30890	138-6529	Hemostasis
LO35041	LOCK-G-26-235-C-N	LO35040	TBD	Hemostasis
LO35045	LOCK-G-26-235-C-N	LO35040	TBD	Hemostasis
MB23031	MBD-AM-BB-2	MB23030	125-2734	Dilation
MB23071	MBD-BM-BB-2	MB23070	125-2736	Dilation
MB23111	MBD-CM-BB-2	MB23110	125-2738	Dilation
MB23151	MBD-DM-BB-2	MB23150	125-2740	Dilation
MB23191	MBD-EM-BB-2	MB23190	125-2742	Dilation
MB23221	MBD-FM-BB-2	MB23220	125-2744	Dilation
ML32651	MBLS-6F-NL	ML32650	138-5348	Hemostasis
NK62721	ENK10	NK62720	139-2221	Polypectomy
NT54511	NT-E-30/23-P	NT54510	133-5543	Polypectomy

Order Number	Ref #	UPN Pouch	Henry Schein Item Number	Product Family
NT54515	NT-E-30/23-P	NT54514	135-6057	Polypectomy
PC54011	BBP-50	PC54010	129-0676	GI Accessories
PC54021	BBP-50-S	PC54020	131-2219	GI Accessories
PC54031	BBP-50-M	PC54030	131-2214	GI Accessories
PS51011	NPFS01-01023230	PS51010	129-0186	Polypectomy
PS51021	NPFS01-01523230	PS51020	129-0673	Polypectomy
PS51031	NPFS01-02023230	PS51030	129-0674	Polypectomy
PS51041	NPFS01-02423230	PS51040	129-0675	Polypectomy
PS51051	NPFS01-03023230	PS51050	129-0187	Polypectomy
PS51071	NPFS02-01023230	PS51070	129-0188	Polypectomy
PS51081	NPFS02-01523230	PS51080	129-0189	Polypectomy
PS51091	NPFS02-02023230	PS51090	129-0190	Polypectomy
PS51101	NPFS02-02423230	PS51100	129-0191	Polypectomy
PS51111	NPFS02-03023230	PS51110	129-0192	Polypectomy
PS51121	NPFS02-03623230	PS51120	129-0193	Polypectomy
PS52031	NPFS02-11523230	PS52030	129-0194	Polypectomy
PS52041	NPFS02-12523230	PS52040	129-0195	Polypectomy
QL22011	MDBR1A0504A	QL22010	131-2220	Infection Prevention
QL22021	MDBR1B0016A	QL22020	131-2222	Infection Prevention
QL22031	MDBR1B0015A	QL22030	131-2221	Infection Prevention
QL22041	MDBR1B0005A	QL22040	131-2215	Infection Prevention
QL22051	MDBR1B0013A	QL22050	131-2216	Infection Prevention
QL22061	MDBR1B0014A	QL22060	131-2217	Infection Prevention
QL22511	MDKT1A0001A	QL22510	140-0706	Infection Prevention
QL22521	MDKT1A0002A	QL22520	140-0705	Infection Prevention
QL22531	MDKT1A0003A	QL22530	140-0557	Infection Prevention
QL22541	MDKT1A0004A	QL22540	140-0704	Infection Prevention
QL22711	MDJL1A0001A	QL22710	140-1637	Infection Prevention
QL22811	MDSB1B0001A	QL22810	140-1682	Infection Prevention
QL26021	MDBB1AY002A	QL26020	TBD	GI Accessories
QL30011	MDHL100000A	QL30010	138-8285	Hemostasis
QL54012	MDPT1A0002A	QL54010	136-8983	GI Accessories
QL70021	MD-A-BV02	QL70020	138-6532	Infection Prevention
QL70051	MD-B-BV02	QL70050	138-6534	Infection Prevention
RC30381	ROCC-F-26-235-C-R	RC30380	128-5655	Hemostasis
RC30385	ROCC-F-26-235-C-R	RC30384	128-5657	Hemostasis
RC30441	ROCC-D-26-235-C-N	RC30440	132-5723	Hemostasis
RC30445	ROCC-D-26-235-C-N	RC30444	132-5724	Hemostasis
RC30541	ROCC-H-26-235-C-N	RC30540	136-6310	Hemostasis
RC30545	ROCC-H-26-235-C-N	RC30544	136-6309	Hemostasis
SC62241	WP-18/2200	SC62240	133-5544	Polypectomy
SC62301	WP-23/2500	SC62300	133-5545	Polypectomy
TP54411	BBPA00001	TP54410	TBD	GI Accessories
TP54421	BBPA10001	TP54420	TBD	GI Accessories



**[MTEndoscopy.com](http://MTEndoscopy.com)**

©2021 Micro-Tech Endoscopy USA Inc. All rights reserved. All other trademarks are the property of their respective owners. MPV001S Rev 5